HQ:CSC38/Cr-Au

AFM Probe with 3 Different Gold Coated Contact Mode AFM Cantilevers


AFM probes of the HQ:CSC38 series have three different soft contact mode AFM cantilevers on one side of the holder chip. They can be used in various applications.

The HQ AFM probes offer high consistency of the AFM tip radius, the AFM cantilever reflectivity and the quality factor.

The overall 30 nm Au coating with 20 nm Cr sublayer is electrically conductive and chemically inert. It also ehnances the laser reflectivity of the AFM cantilevers in air and liquids. The resulting coated AFM tip radius is below 35 nm. The coating may cause AFM cantilever bending up to 3°.

Coating

Gold Overall

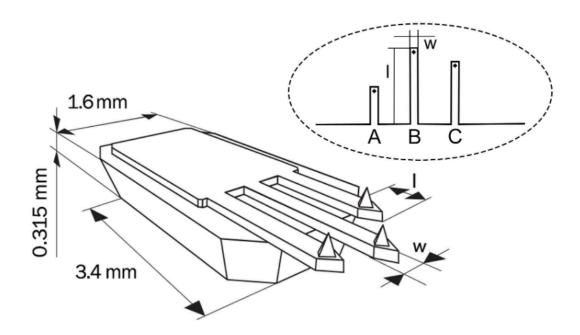

AFM Probe Specifications

AFM Tip

SHAPE	HEIGHT	FULL CONE ANGLE	RADIUS
Rotated	15 μm (12 – 18 μm)*	40°	< 35 nm

AFM Cantilever

CANTILEVER	SHAPE	FORCE CONST.	RES. FREQ.	LENGTH	WIDTH	THICKNESS
Cantilever A	Beam	0.09 N/m (0.01 – 0.36 N/m)*	20 kHz (8 – 32 kHz)*	250 μm (1 – 255μm)*	32.5 μm (29.5 – 35.5μm)*	1μm (0.5 – 1.5 μm)*
Cantilever B	Beam	0.03 N/m (0.003 – 0.13 N/m)*	10 kHz (5 – 17 kHz)*	350 μm (1 – 355μm)*	32.5 μm (29.5 – 35.5μm)*	1μm (0.5 – 1.5 μm)*
Cantilever C	Beam	0.05 N/m (0.005 – 0.21 N/m)*	14 kHz (6 – 23 kHz)*	300 μm (1 – 305μm)*	32.5 μm (29.5 – 35.5μm)*	1μm (0.5 – 1.5 μm)*

^{*} typical values